

## GP-30T

### Pneumatic grip 30kN

GP-30T are supplied in pairs, but can also be ordered separately.


<b>Item No.</b>	<b>GP-30T</b>
Clamping force	30 kN/7 bar
Temperature range	0°C to +70°C Other temperature ranges on request
Body	Aluminum, anodized
Pneumatic connection	1/8"
Weight per grip	10kg each grip (without jaws)
Scope of delivery	1 pair of grips (without jaws)



**Jaws for GP-30T:** Scope of delivery 1 set = 4 jaws

Item No.	Jaw surface	Clamping surface H x W	Opening	Weight / pair
GP-30T-B	Blank (smooth) jaws	40x120mm	0-12 mm *	1.08kg
GP-30T-BG	Rubber jaws	40x120mm	0-10 mm *	1.09kg
GP-30T-BP	Pyramid (serrated) jaws	40x120mm	0-12 mm *	1.04kg
GP-30T-BV6	V-jaws	40x120mm	6-17 mm Ø *	0.83kg
GP-30T-BW	Wave jaws	40x120mm	0-12 mm *	1.02kg
GP-30T-BD	Diamond jaws	40x120mm	0-12 mm *	1.08kg
GP-30T-BW200	Wave jaws	60x200mm	0-8 mm	3.61kg
GP-30T-BW210	Wave jaws	60x210mm	0-8 mm	3.76kg

\*Opening specs can be changed upon request: 4-16mm, 8-20mm and 10-22 mm (depending on jaw type)  
Please let us know prior to ordering; modification after delivery is not possible.

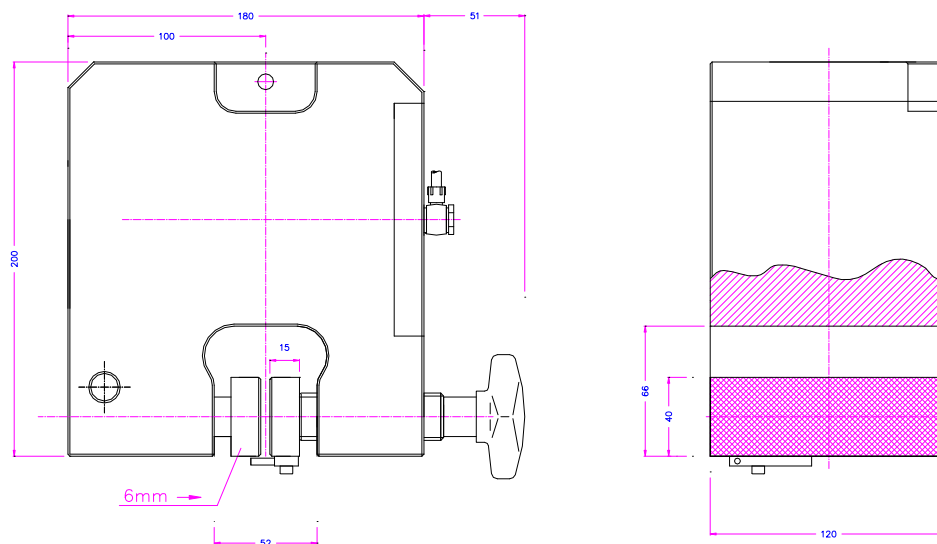

## GP-30T

Other versions and customer-built solutions



### Jaws Surfaces:



### Accessories:



#### GP-1 Non-Locking Foot Switches

2 pneumatic foot switches, tubes and fittings

Non-locking: The jaws of the grip open by pressing the foot pedal and close after removing the foot from the pedal.

Scope of delivery: 1 pair of foot switches including tubes and fittings



#### GP-2 Locking Foot Switches

2 pneumatic foot switches, tubes and fittings,

With locking function: The jaws of the grip open by pressing the foot pedal and remain opened after removing the foot from the pedal. The jaws close by pressing the pedal one more time.

Scope of delivery: 1 pair of foot switches including tubes and fittings



#### GP-PR

Optional for GP-1 or GP-2

Pressure regulator

To adjust air pressure: 1-10 bar



#### GP-3 Locking Hand Switches

1 unit of pneumatic hand switches, tubes and fittings

With locking function: Handle with 2 positions. The jaws of the grip close by pulling the handle and remain closed. The jaws open by pushing the handle back.

Scope of delivery: 1 unit of hand switches including tubes and fittings



#### GP-4 Locking Hand switches

The hand switches are directly mounted on the grips. 3/2 directional valve with locking function.

Scope of supply: 2 hand switches including tubes and fittings



#### GP-CR Compressor for laboratories, oil-lubricated

Max pressure 16 bar, 230V, noise level: 38dB(A)/1m

Displacement 5 l/min, tank size 3.1 litre (0.82 gallon)

Dimensions: 340x290x340 mm (L x W x H); weight 16 kg

Different voltages available on request